

RAUE
Kosmetik- & Fußpflegebedarf

Datum
09-04-19

Start
09:28:48


Stopp
09:47:33

S/N
ST01-00260

Programm
134

Prozessnummer
00009

Firmware
27.3

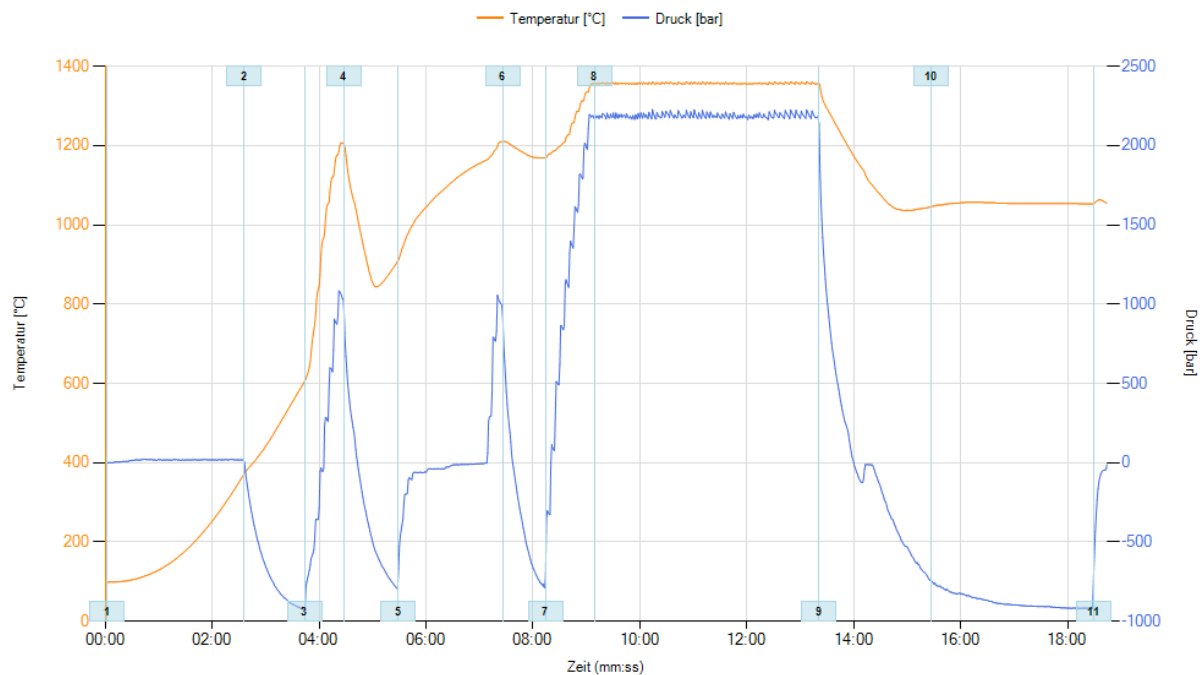

	Phase	Zeit [mm:ss]	Schritt [mm:ss]	p [bar]	t [°C]	t calc. [°C]
0	Locking door	00:01		0.002	9.9	9.9
1	Heating up	00:02	02:33	0.019	36.9	36.9
2	Vacuum I	02:36	01:07	-0.925	60.3	60.3
3	Steam Pulse I	03:44	00:43	1.031	120.5	120.5
4	Vacuum II	04:28	01:00	-0.795	90.8	90.8
5	Steam Pulse II	05:29	01:56	0.996	121.0	121.0
6	Vacuum III	07:26	00:47	-0.789	116.9	116.9
7	Pressurizing	08:14	00:54	2.174	135.7	135.7
8	Sterilizing	09:09	04:11	2.175	135.7	135.6
9	Depressurizing	13:21	02:05	-0.740	104.5	104.5
10	Drying	15:27	03:01	-0.918	105.3	105.3
11	Equalizing	18:29	00:14	-0.043	105.6	105.6
12	Ulocking door	18:44		-0.011	105.5	105.5

Sterilization Complete

Sterilisationstemper 134

Zeit 04:11

Min. Temperatur 135,6

Max. Temperatur 136

Min. Druck 2,16

Max. Druck 2,23

Teststerilisation bei Lieferung
Inhalt: unverpackt
Gelenkzange
Pinzette
Tamponierer

Label